


MetaClip-P & Appliers

Polymer Ligating Clips & Appliers



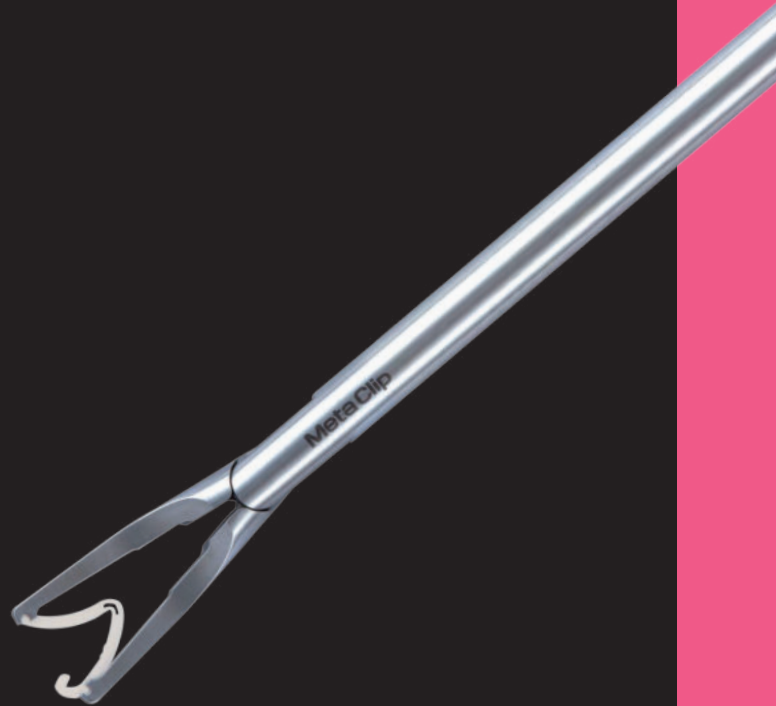

meta Surgical


MetaClip-P & Appliers

Polymer Ligating Clips & Appliers




– MetaClip-P Features

- Non slip type
- Clear Locking sound
- Biocompatible Polymer
- Non-absorbable Polymer
- Flexible Hinge
- High tensile strength

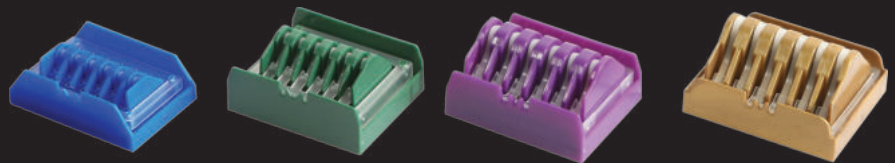

– Clip Specification

Vessels / Tissues bundle size range

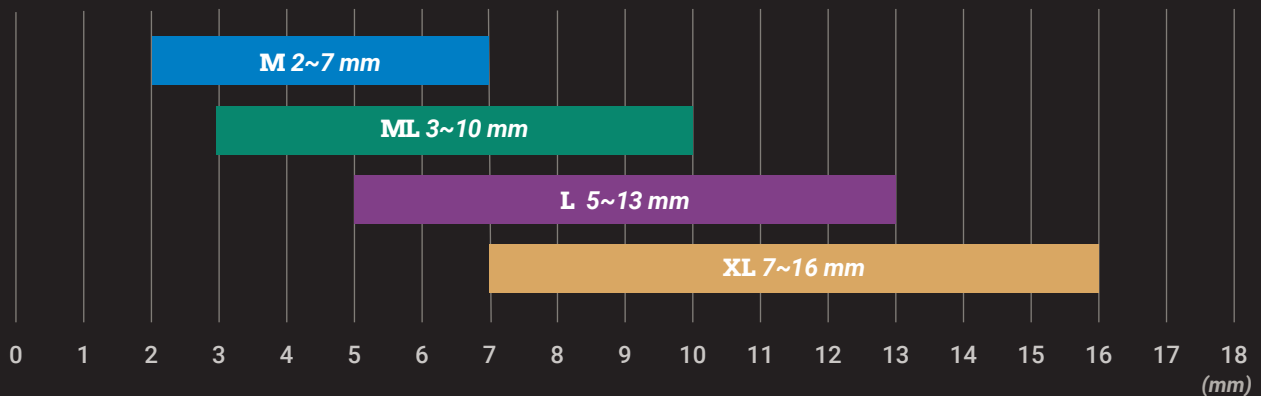


| Clip | Ref No. | Clips | Cartridges / Box | Clips / Box |
|------|---------|-------|------------------|-------------|
| M | PCMM01A | 6 | 14 | 84 |
| ML | PCML02A | 6 | 14 | 84 |
| L | PCLL03A | 6 | 14 | 84 |
| XL | PCXL04A | 6 | 14 | 84 |



CLIP

Non-Absorbable Polymer

is a biocompatible polymer that does not interfere with CT, MRI or X-ray.

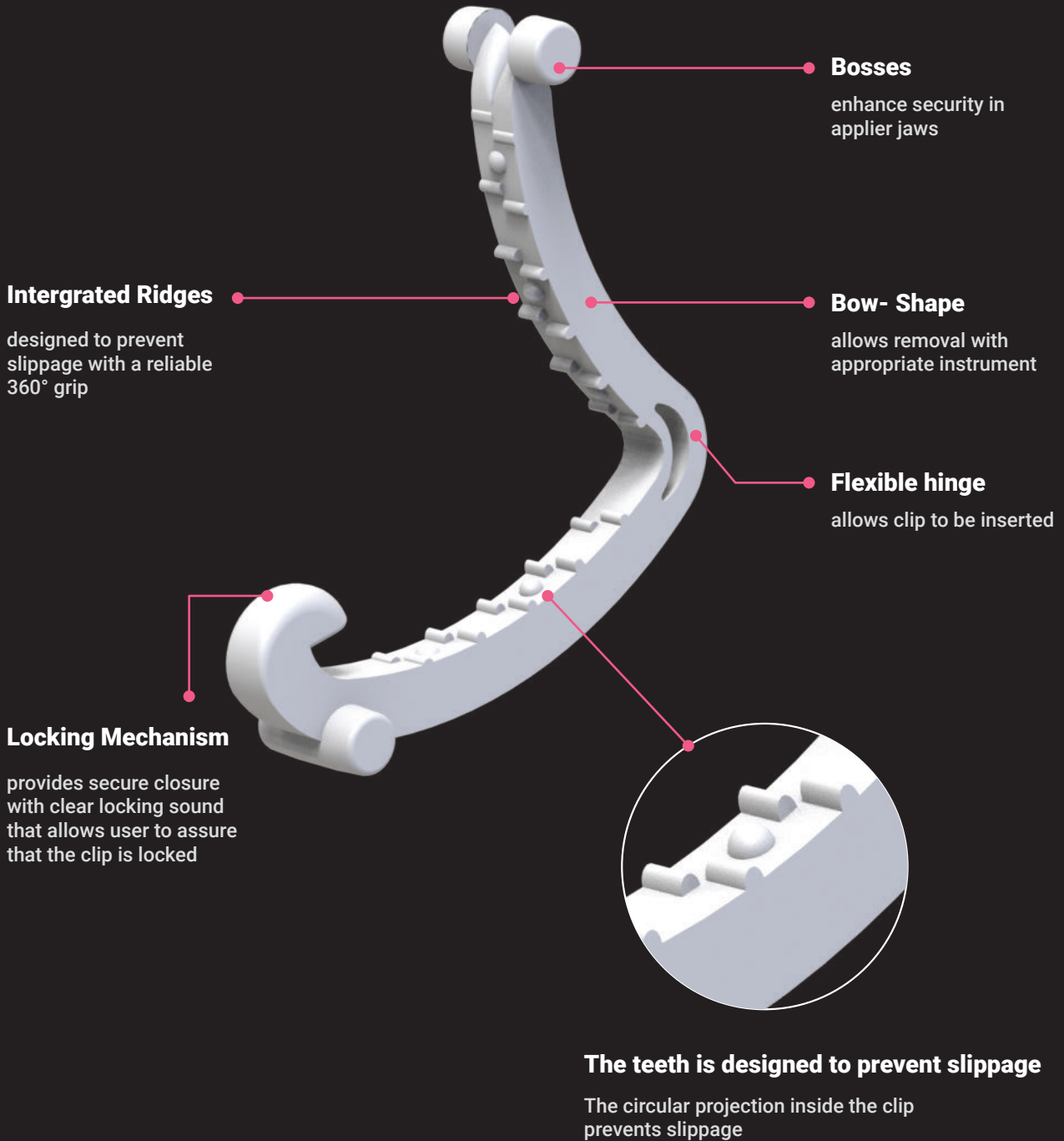



Appliers

Proper Jaw Angle

for better loading it is necessary to position the applier jaw perpendicular to the cartridge body



– Appliers Features

- Non - detachable design for easier maintenance and durability
- Made of high quality stainless steel
- Color coded size matching cartridge
- Compatible with corresponding sizes of Meta Clip-P polymer



– Applier Sterilization

○ Sterilization Conditions

Wrap with a disinfection Cloth.


○ Instrument State

Pre-vacuum

| APPLIER | TEMPERATURE | EXPOSURE TIME | DRYING TIME |
|------------------|-------------|---------------|-------------|
| All Type Applier | 134-135° C | 3 minutes | 30 minutes |




◦ OPEN APPLIER

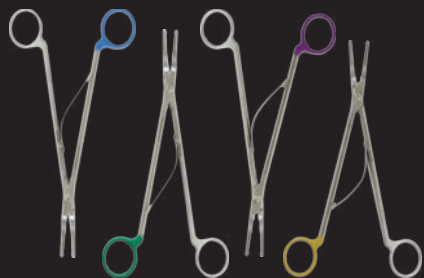

| Clip | Model No | Length |
|------|----------|--------|
| M | OMAPMMA | 280mm |
| ML | OMAPMLA | 280mm |
| L | OMAPLLA | 280mm |
| XL | OMAPXLA | 280mm |

◦ ENDO APPLIER

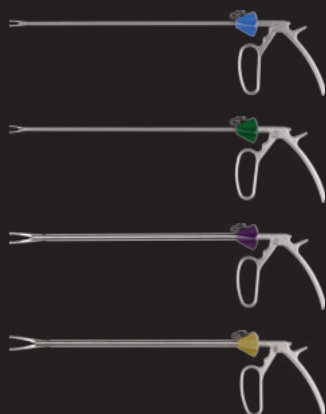


| Clip | Model No | Shaft Length | Shaft Diameter |
|------|----------|--------------|----------------|
| M | EMAPMMA | 330mm | 5mm |
| ML | EMAPMLA | 330mm | 5mm |
| L | EMAPLLA | 330mm | 10mm |
| XL | EMAPXLA | 330mm | 10mm |

◦ OPEN REMOVER



| Clip | Model No | Length |
|------|----------|--------|
| M/ML | ORAPMLA | 280mm |
| L/XL | ORAPXLA | 280mm |

◦ ENDO REMOVER

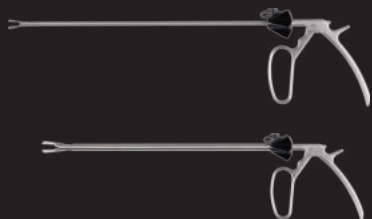

| Clip | Model No | Shaft Length | Shaft Diameter |
|------|----------|--------------|----------------|
| M/ML | ERAPMLA | 320mm | 5mm |
| L/XL | ERAPXLA | 320mm | 10mm |



◦ Order Information

| Clip | Open Applier | Endo Applier | Open Remover | Endo Remover |
|------|--------------|--------------|--------------|--------------|
| M | OMAPMMA | EMAPMMA | ORAPMLA | ERAPMLA |
| ML | OMAPMLA | EMAPMLA | | |
| L | OMAPLLA | EMAPLLA | ORAPXLA | ERAPXLA |
| XL | OMAPXLA | EMAPXLA | | |



meta Surgical

Head Office

Meta Biomed Co., Ltd.

270 Osongsaengmyeong1-ro,
Osong-eup, Heungdeok-gu,
Cheongju-si, Chungbuk, Korea

Phone: +82 43 230 8846

Fax: + 82 43 230 8969

boram@meta-biomed.com

www.meta-biomed.com

Meta Biomed Inc.

121 Independence Lane
Chalfont, PA 18914, USA

Phone: +1 267 282 5893

america@metabiomed-inc.com

www.buymetabiomedonline.com

Meta Biomed Europe

Haus der Wirtschaft, Raum 3.11
Wiesenstraße 35
45473 Mülheim an der Ruhr, Germany

Phone: +49 208 309 9190

europe@metabiomed-inc.com

www.buymetabiomedonline.com